



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA
ELECTRONIC ENGINEERING SRL**

Patient No. 1

The picture "1 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...
during the period between ... 10/01/2005..... and 10/02/2005.....

Electrodes under the following configurations have been used:

- Configuration 1A (electrodes A, B, C, D)

for each electrode it has been set: density of current ...2.4..... mA / cmq
 density of power ... 100.....mWatt / cmq

- Configuration 2C (electrodes A, B, C, D)

for each electrode it has been set: density of current ...2.4.....mA / cmq
 density of power... 100.....mWatt / cmq


For each configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "1 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.


Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine


.....

Files enclosed:

- Pictures "1 Before" and "1 After"
- Drawings of the electrodes configurations as used


TELEA ELECTRONIC ENGINEERING S.r.l.
Via Leopardi, 23 - 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa - VI-
Uff. Commerciale: +390444239519 - +390444918179 - Uff. Amministrazione: +390444919322
FAX: +390444914319 E-mail: info@vesalius.it - <http://www.vesalius.it>
Codice Fiscale e Partita I.V.A. 01969870243 - Cap. Soc. Euro 63.191,08 I.V. - Reg. Impr. N. 21604

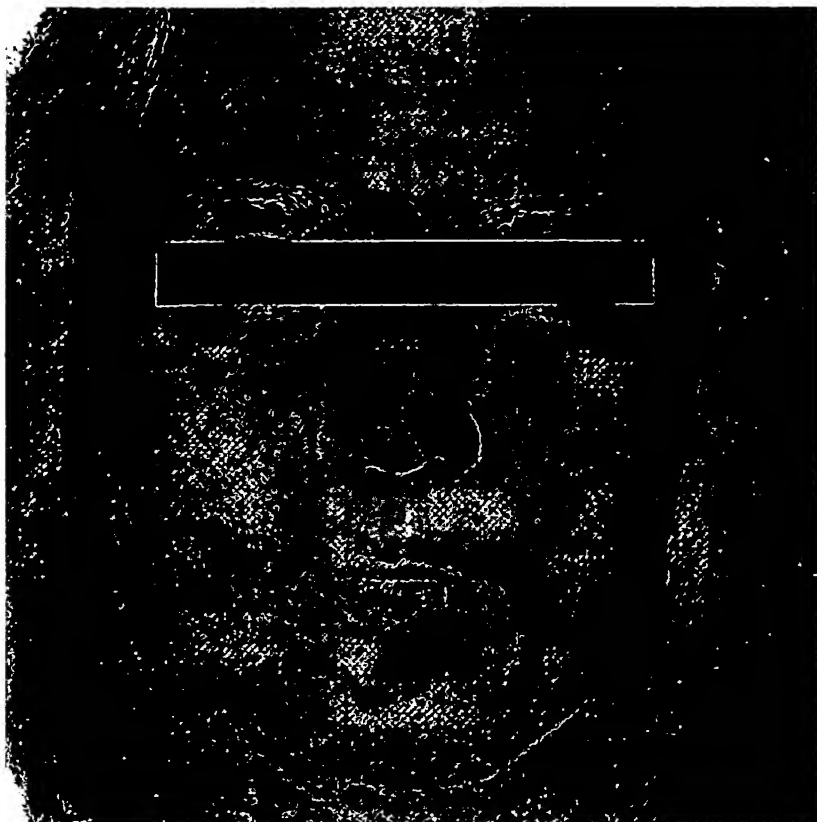

flexon age


1

Before

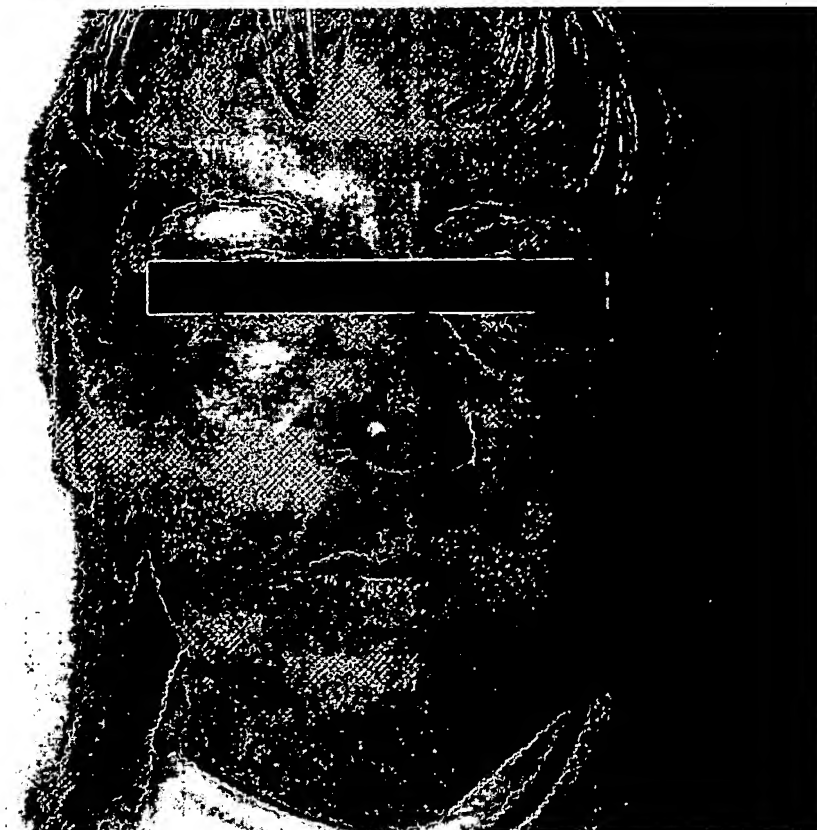


After


Electrodes setting schemes
ELECTRODE nr. 1 (Eye Contour)

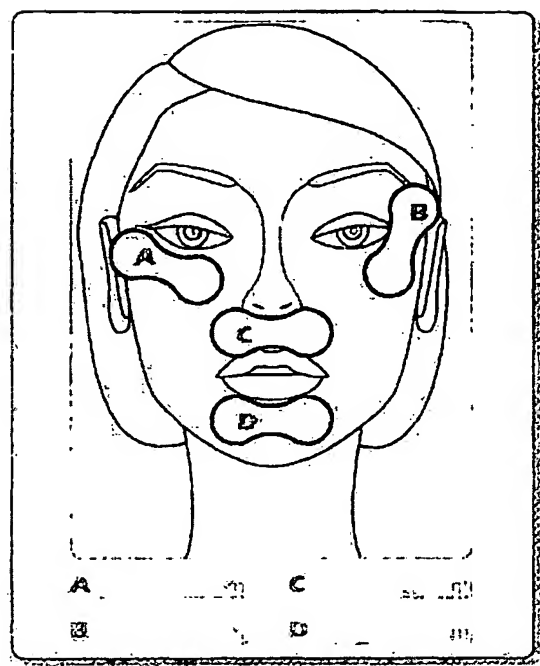

1A

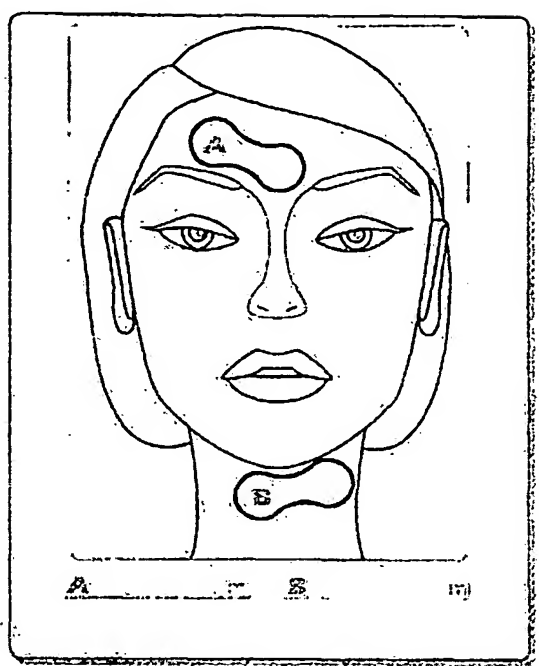

1B

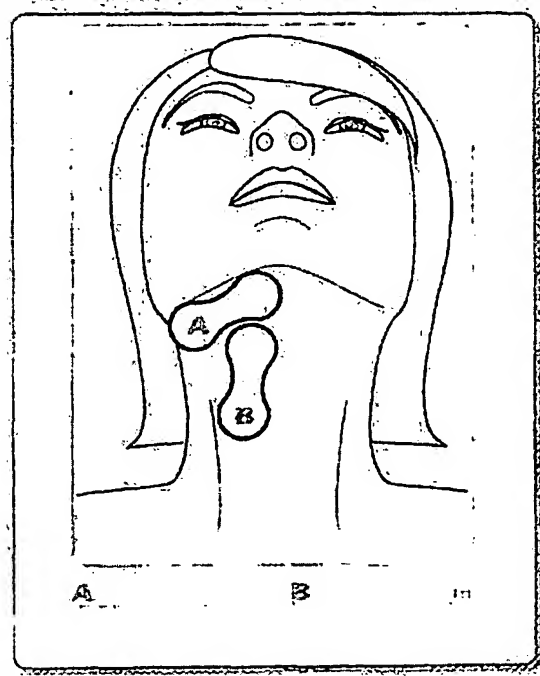


1C


Number of cycles per session: _____

Application time per each electrode: _____


Electrodes setting schemes
 ELECTRODE nr. 2 (Medium Size)

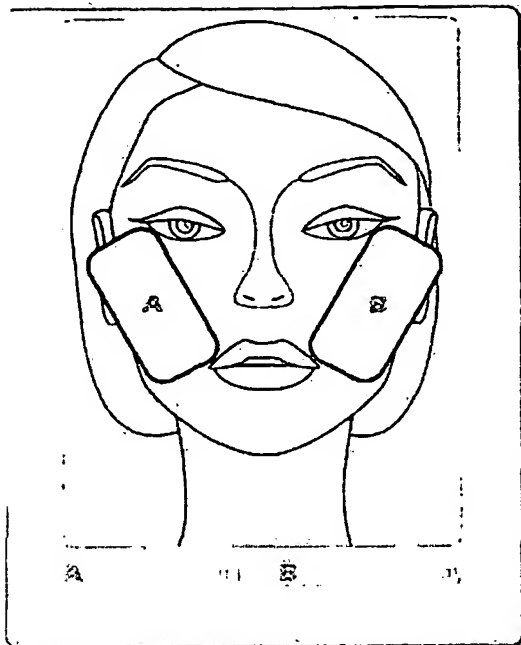


2A

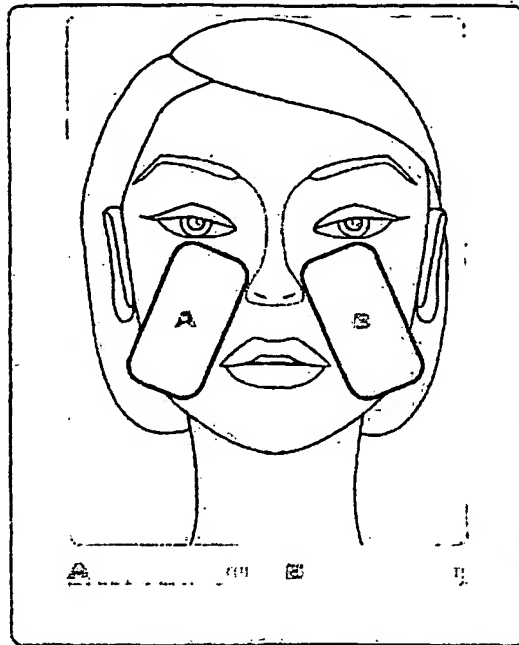

2B

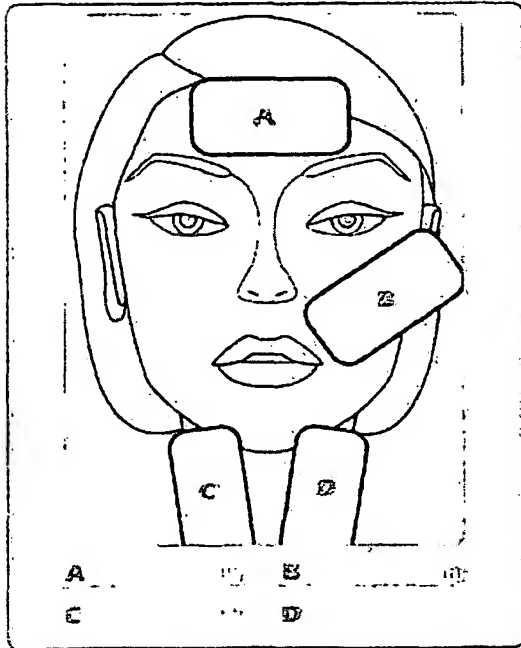

2C

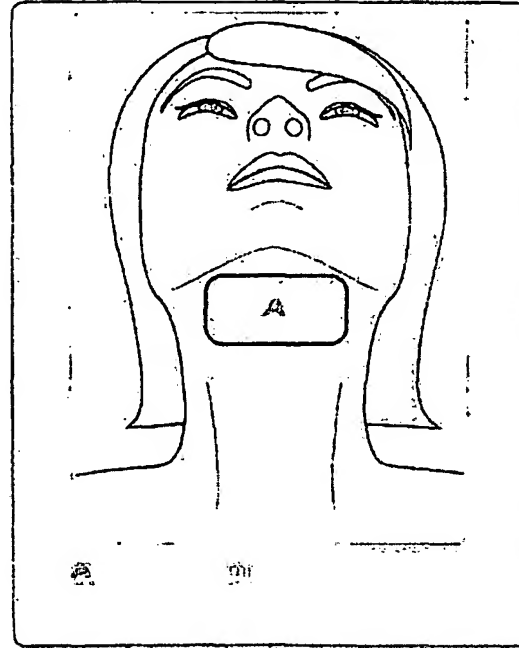

Number of cycles per session: _____

Application time per each electrode: _____



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA
ELECTRONIC ENGINEERING SRL**

Patient No. 2

The picture "2 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...
during the period between ...10/01/2005..... and10/02/2005.....

Electrodes under the following configuration have been used:

- Configuration 1A (electrodes A, B, C, D)

for each electrode it has been set: density of current ...2.4..... mA / cmq
 density of power ...100.....mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The pictures "2 After" have been taken at the end of the therapy, and they show considerable improvement with respect to the initial pictures.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine



.....

Files enclosed:

- Pictures "2 Before" and "2 After"
- Drawings of the electrodes configurations as used


REGION 8

2


Before

After


Electrodes setting schemes
ELECTRODE nr. 1 (Eye Contour)

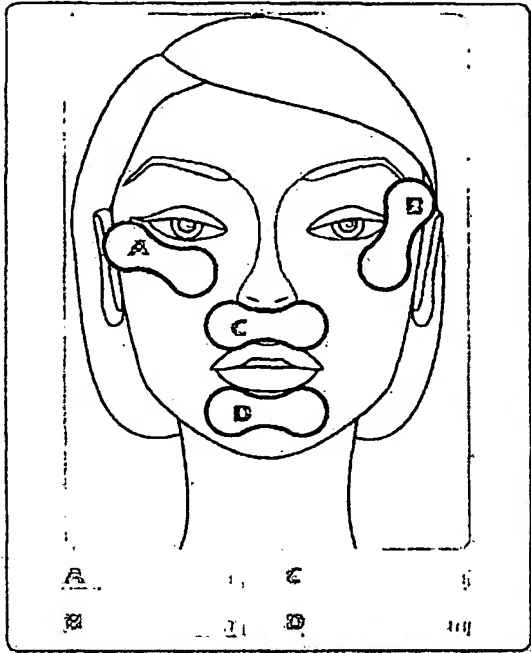

1A

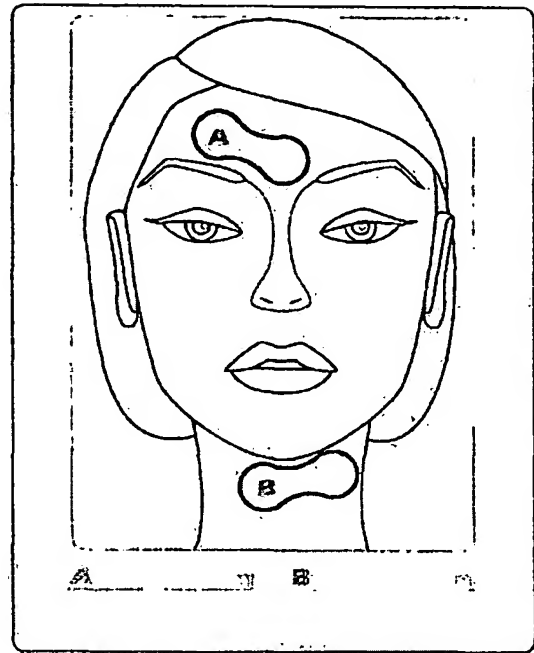

1B

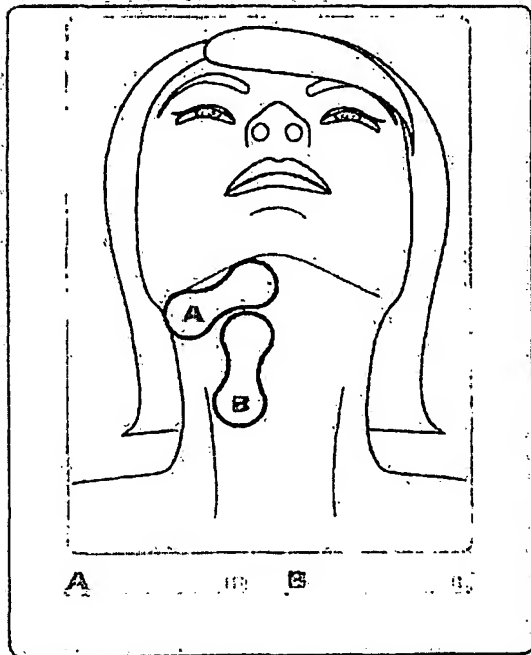

1C

Number of cycles per session: _____

Application time per each electrode: _____



ELETRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

Patient No. 3

The picture "3 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...
during the period between ...10/01/2005... and ...10/02/2005...

Electrodes under the following configuration have been used:

- Configuration 1A (electrodes A, B)

for each electrode it has been set: density of current ...2.4... mA / cmq
density of power ...100... mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.


The picture "3 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine

Files enclosed:

- Pictures "3 Before" and "3 After"
- Drawings of the electrodes configurations as used

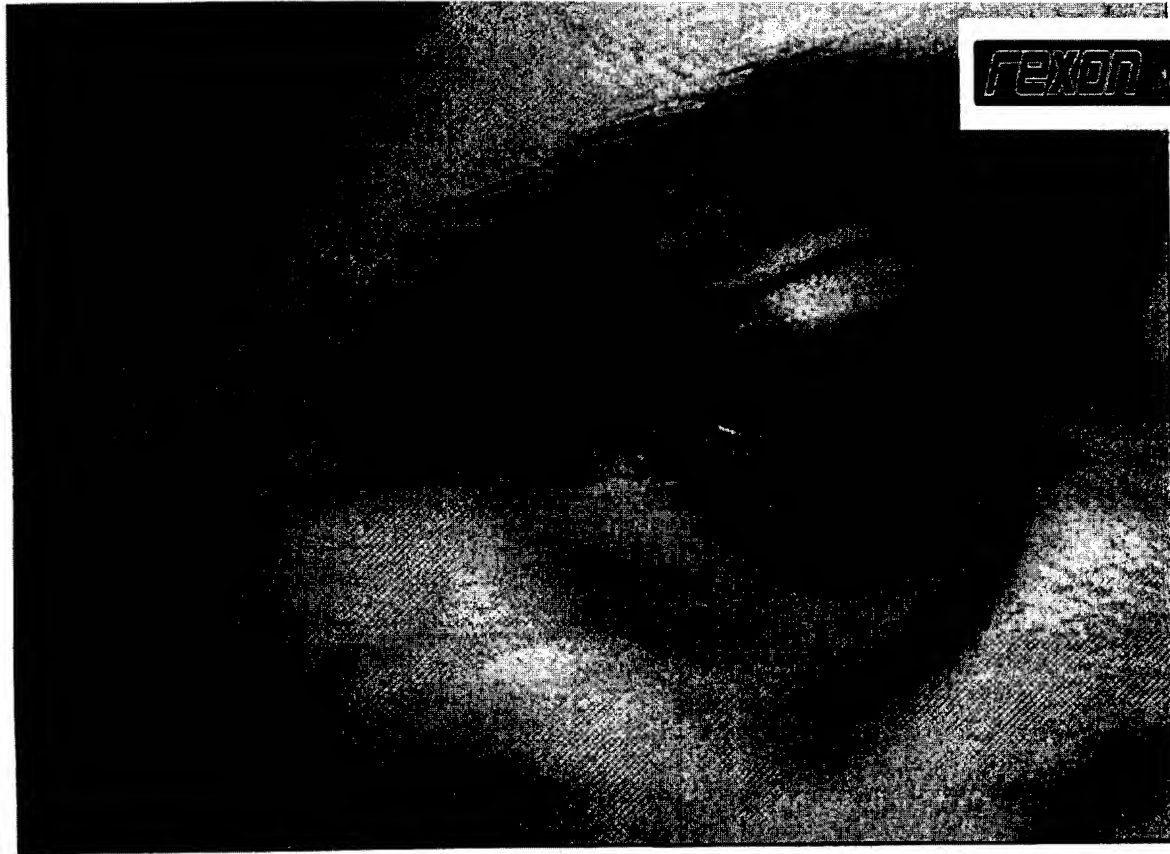


3

Before




After


Electrodes setting schemes
ELECTRODE nr. 1 (Eye Contour)

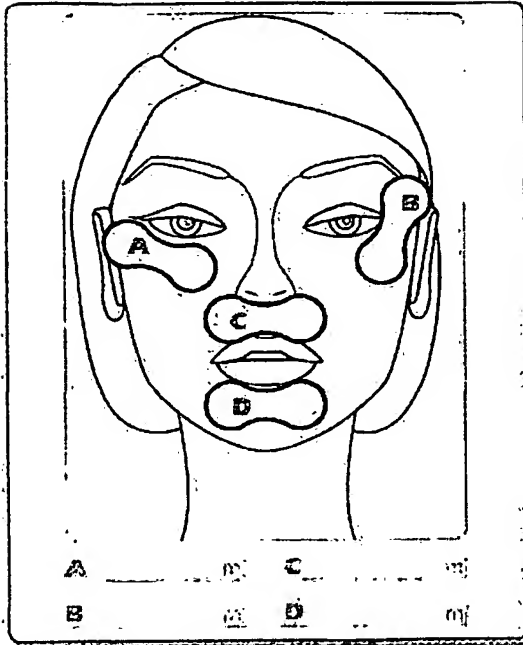

1A

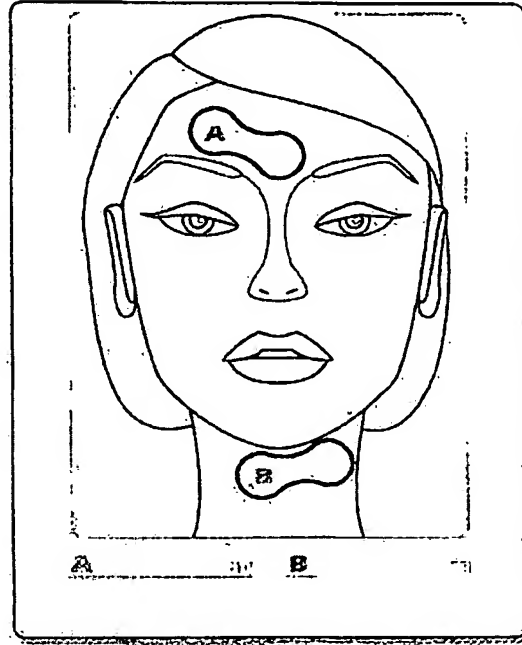

1B

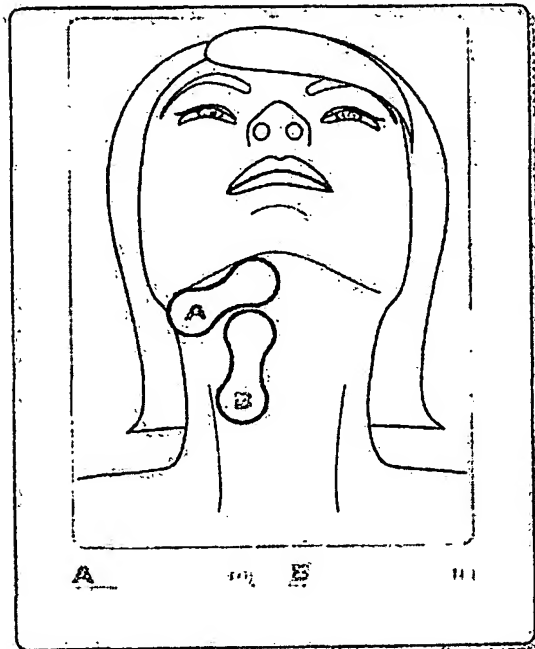

1C

Number of cycles per session:

Application time per each electrode:



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

Patient No. 4

The picture "4 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106... during the period between ... 10/01/2005..... and10/02/2005.....

Electrodes under the following configuration have been used:

- Configuration 3A (electrode C)

the electrode has been set: density of current ...2.4..... mA / cmq
density of power ... 100.....mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "4 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine.

Files enclosed:


- Pictures "4 Before" and "4 After"
- Drawings of the electrodes configuration as used

Page 4

4

Before

After

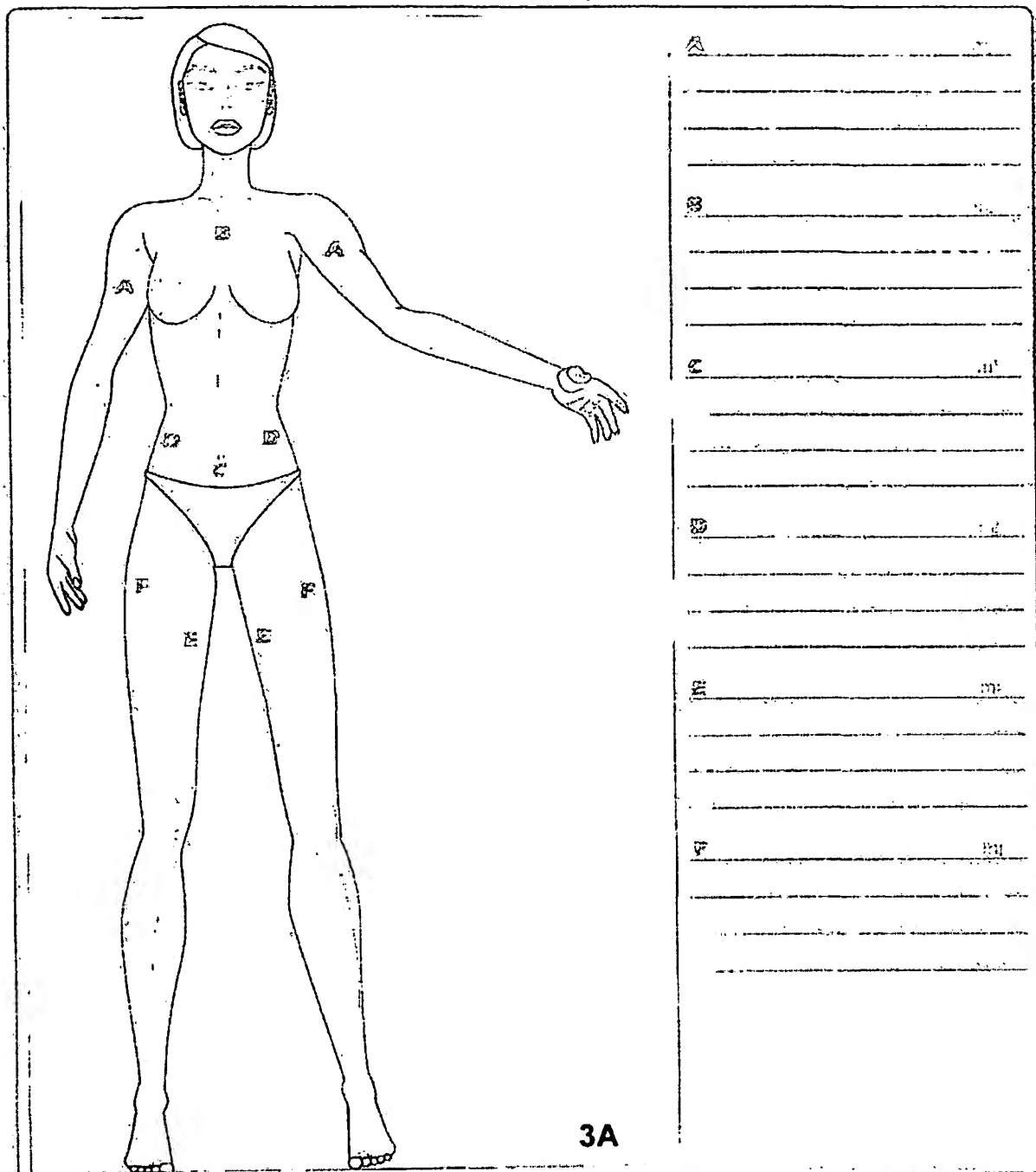

Number of cycles per session: _____

Application time per each electrode: _____



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

Patient No. 5

The picture "5 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106... during the period between ...15/02/2005..... and15/03/2005.....

Electrodes under the following configuration have been used:

- Configuration 3A (electrode C)

the electrode has been set: density of current ...2.4..... mA / cmq
density of power ...100..... mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "5 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.


The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine

Files enclosed:

- Pictures "5 Before" and "5 After"
- Drawings of the electrodes configuration as used

TELEA ELECTRONIC ENGINEERING S.r.l.
Via Leopardi, 23 - 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa - VI-
Uff. Commerciale: +390444239519 - +390444918179 - Uff. Amministrazione: +390444918322
FAX: +390444914319 E-mail: info@vesalius.it - <http://www.vesalius.it>
Codice Fiscale e Partita I.V.A. 01969870243 - Cap. Soc. Euro 63.191,08 I.V. - Reg. Impr. N. 21604

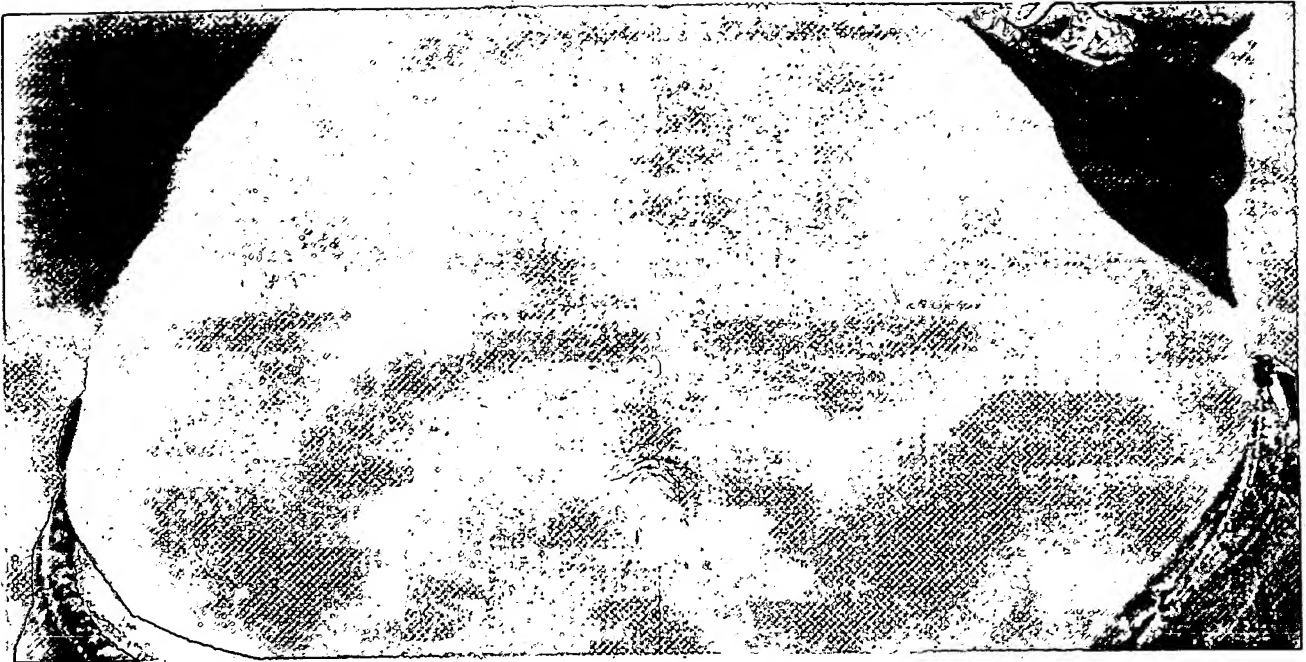

Before



After



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

Patient No. 6

The picture "6 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106... during the period between ...15/02/2005... and ...15/03/2005...

Electrodes under the following configuration have been used:

- Configuration 3A (electrode C)

the electrode has been set: density of current ...2.4... mA / cmq
density of power ...100... mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "6 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine

Files enclosed:

- Pictures "6 Before" and "6 After"
- Drawings of the electrodes configuration as used


REXON AGE

6

Before

After

Electrodes setting schemes
ELECTRODE nr. 3 (Body-front)

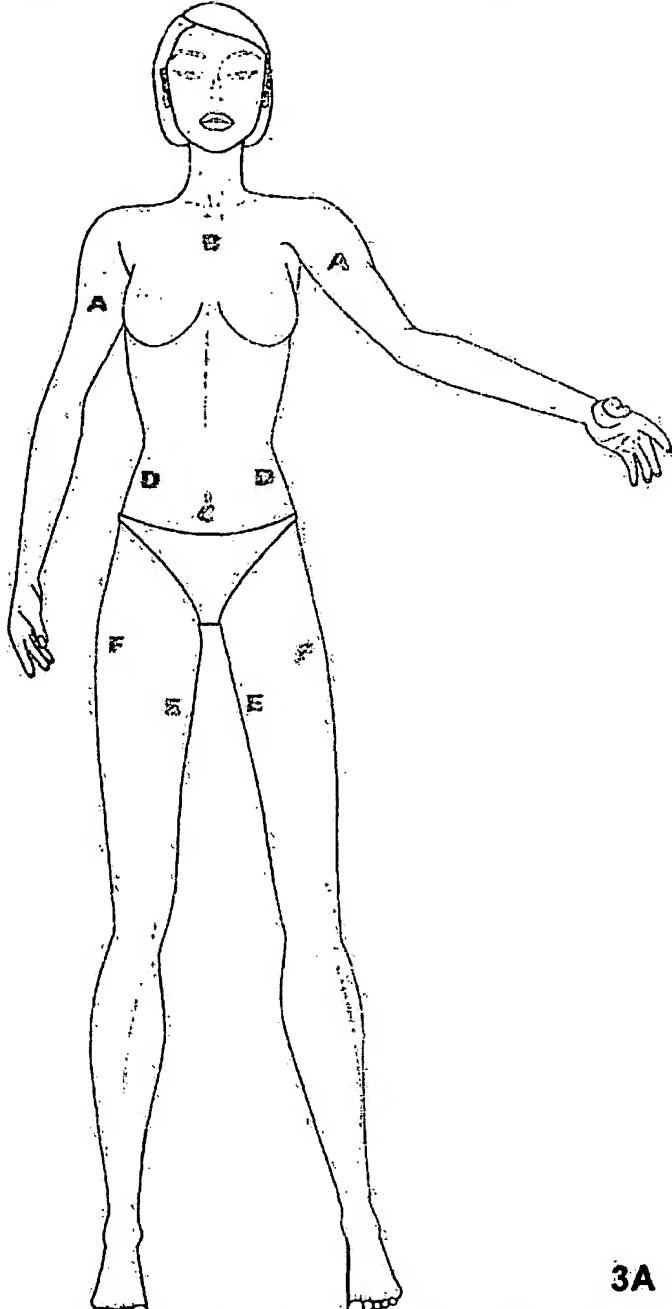

3A

A	_____

B	_____

C	_____

D	_____

E	_____

F	_____

Number of cycles per session: _____
Application time per each electrode: _____



ELETTRONICA
BIOMEDICA
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 45001

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA
ELECTRONIC ENGINEERING SRL**

Patient No. 7

The picture "7 Before" shows a section of a portion of skin and lower tissue removed from the inner part of the left arm of the patient No. 7, to whom the electrode A, according to the configuration 3 A, has then been used.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...
during the period between ...15/02/2005... and ...15/03/2005...

For the electrode A it has been set: density of current ...2.4... mA / cmq
density of power ...100... mWatt / cmq

The electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "7 After" shows a section of a portion of skin and lower tissue removed from the inner part of the left arm of patient No. 7 at the end of the therapy.


From the picture "7 After" we can see that there has been an important improvement in the quantity of collagen and elastine, due to an increased stimulation on the fibroblasts from the energy transmitted.
As a result of the therapy the patient's skin obtained a consistent tightening up and a visible reduction of wrinkles.

Dr. Claudio Trainotti
Doctor and Lecturer in Aesthetic Medicine

Files enclosed:

- Pictures "7 Before" and "7 After"
- Drawings of the electrodes configuration as used


TELEA ELECTRONIC ENGINEERING S.r.l.
Via Leopardi, 23 - 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa - VI -
Uff. Commerciale: +390444239519 - +390444918179 - Uff. Amministrazione: +390444919322
FAX: +390444914319 E-mail: info@vesalius.it - <http://www.vesalius.it>
Codice Fiscale e Partita I.V.A. 01969870243 - Cap. Soc. Euro 63.191,08 I.V. - Reg. Impr. N. 21604

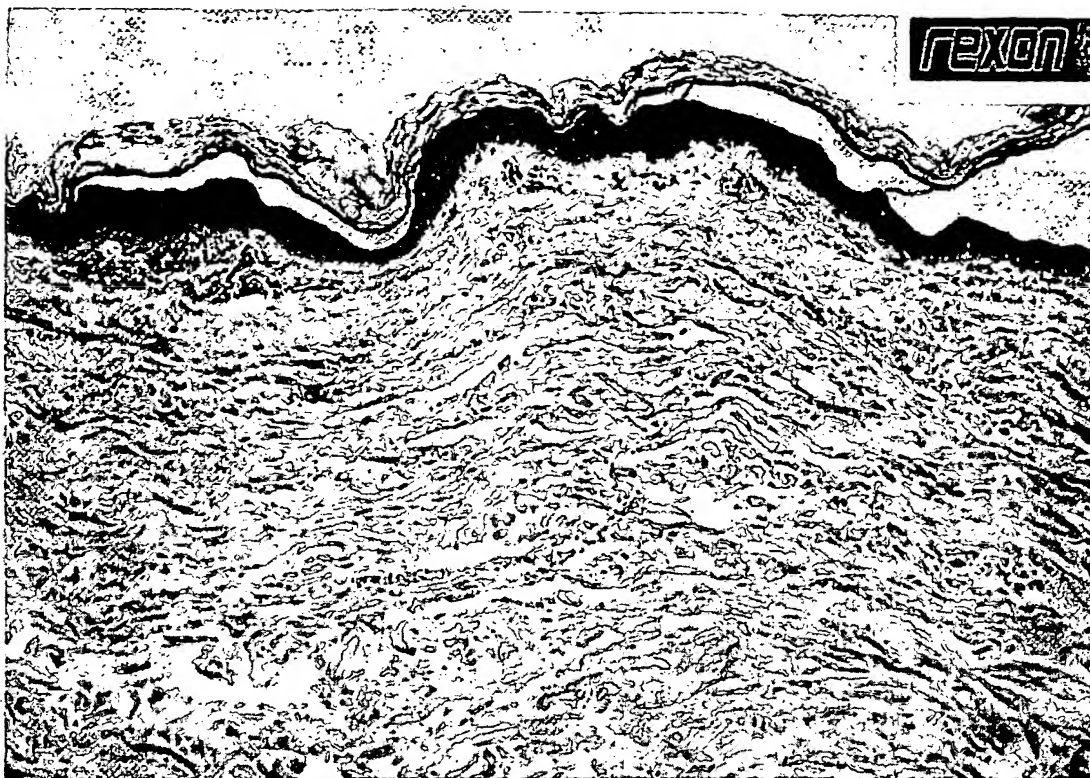

REXON AGE

7


Before

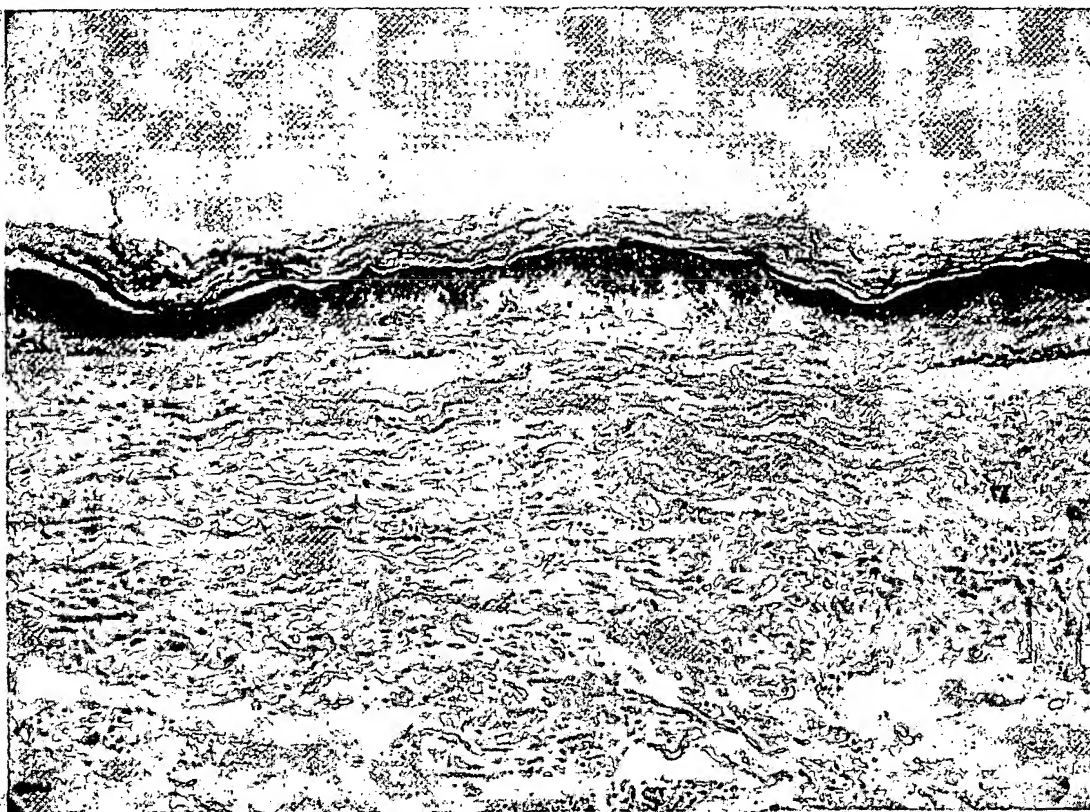

After

Electrodes setting schemes

ELECTRODE nr. 3 (Body-front)

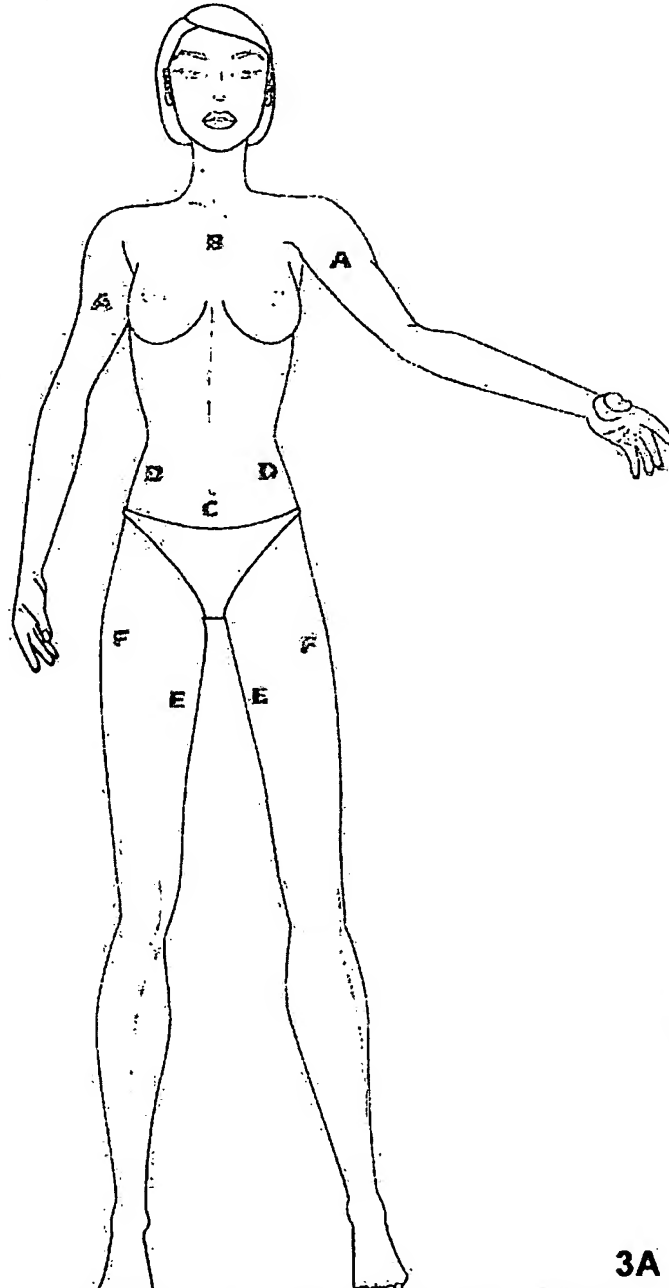

A. _____

B. _____

C. _____

D. _____

E. _____

F. _____

3A

Number of cycles per session: _____

Application time per each electrode: _____