



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

## Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

### Patient No. 1

The picture "1 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...  
during the period between ... 10/01/2005..... and ..... 10/02/2005.....

Electrodes under the following configurations have been used:

- Configuration 1A (electrodes A, B, C, D)

for each electrode it has been set: density of current ... 2.4 ..... mA / cmq  
density of power ... 100 ..... mWatt / cmq

- Configuration 2C (electrodes A, B, C, D)

for each electrode it has been set: density of current ... 2.4 ..... mA / cmq  
density of power ... 100 ..... mWatt / cmq

For each configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.


The picture "1 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

**Files enclosed:**


- Pictures "1 Before" and "1 After"
- Drawings of the electrodes configurations as used

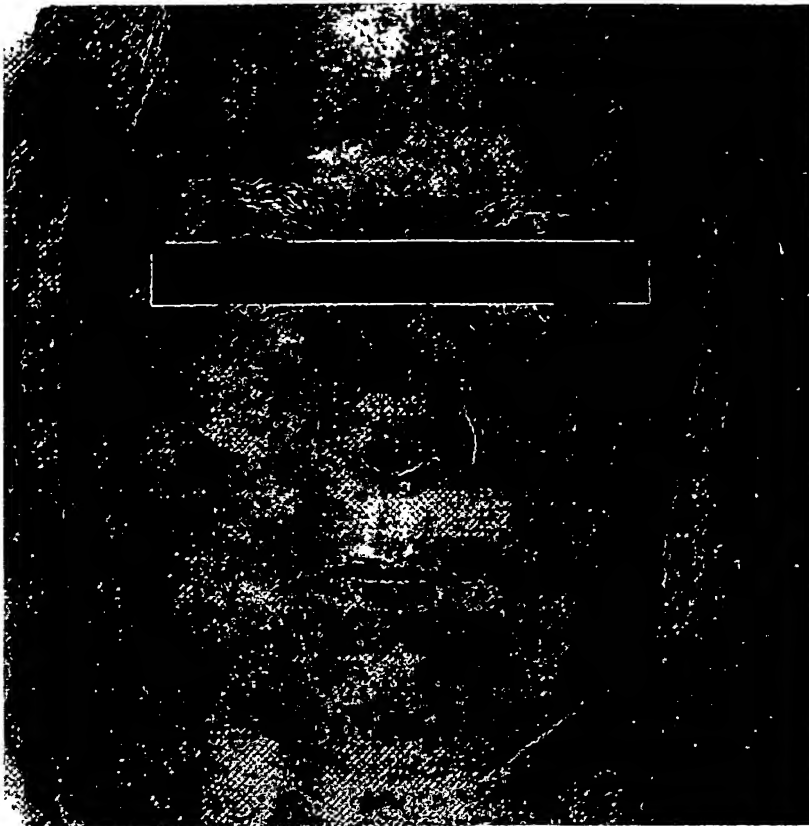

**REXON AGE**




**1**

**Before**

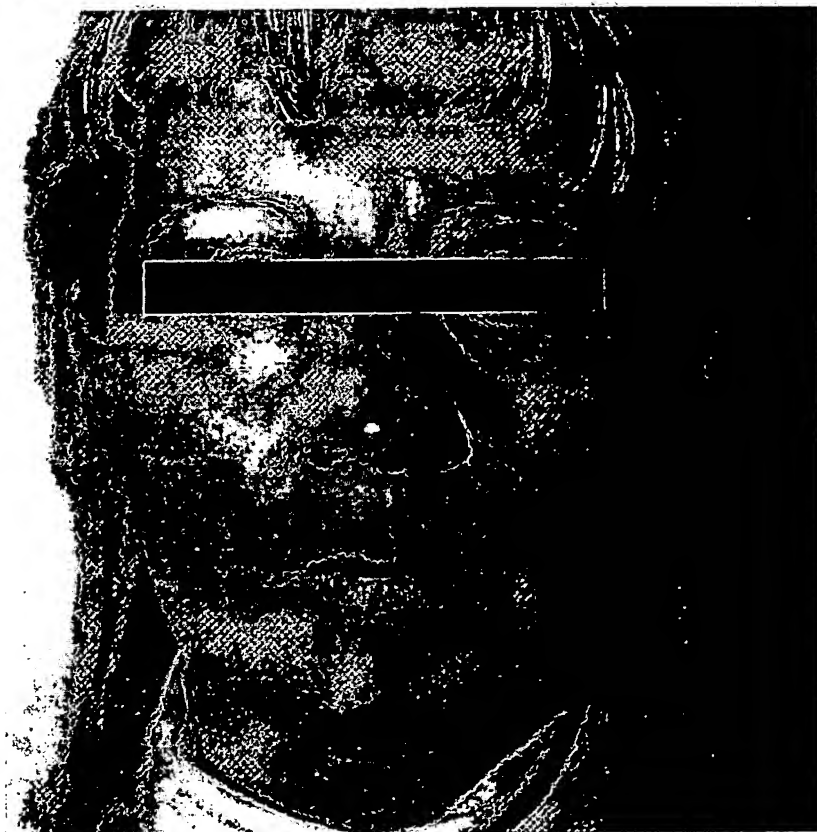

**After**


# Electrodes setting schemes ELECTRODE nr. 1 (Eye Contour)

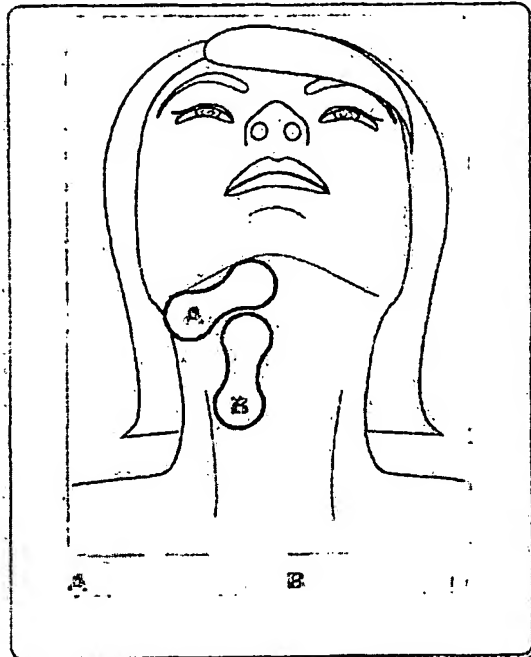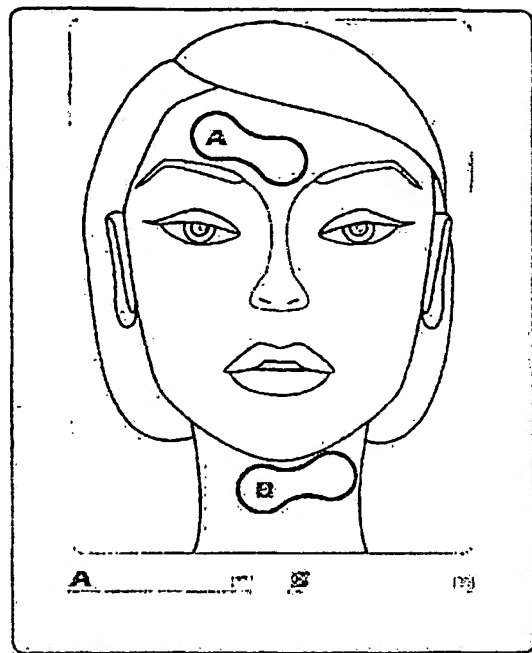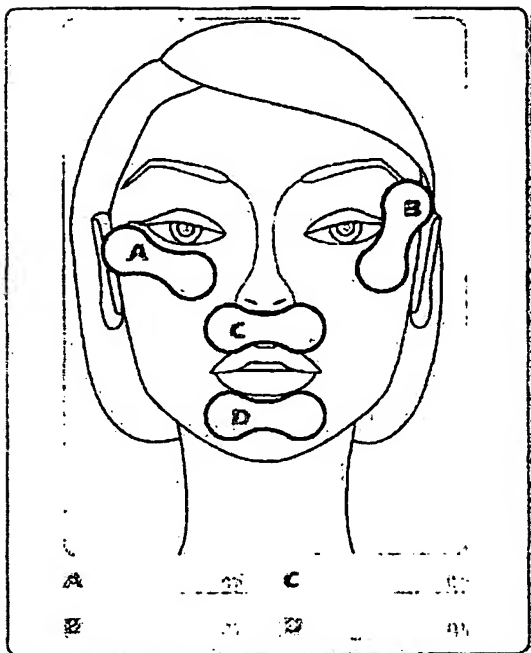


Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_


Electrodes setting schemes  
ELECTRODE nr. 2 (Medium Size)

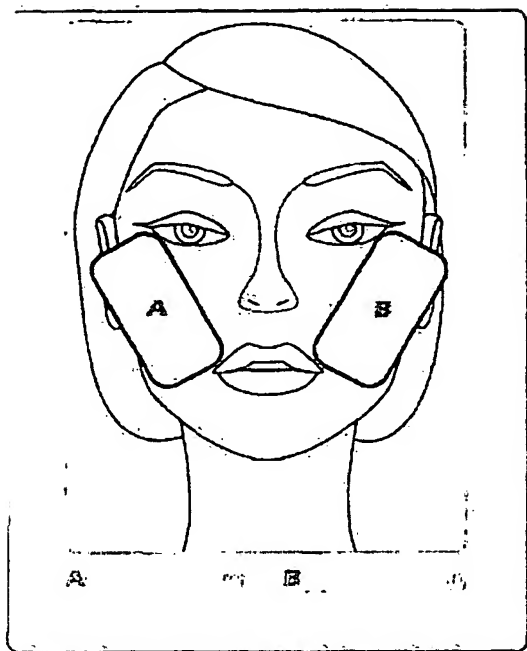


2A

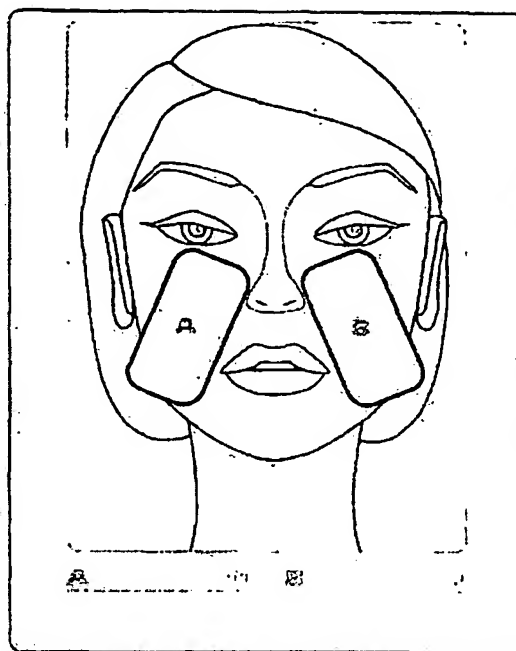

2B

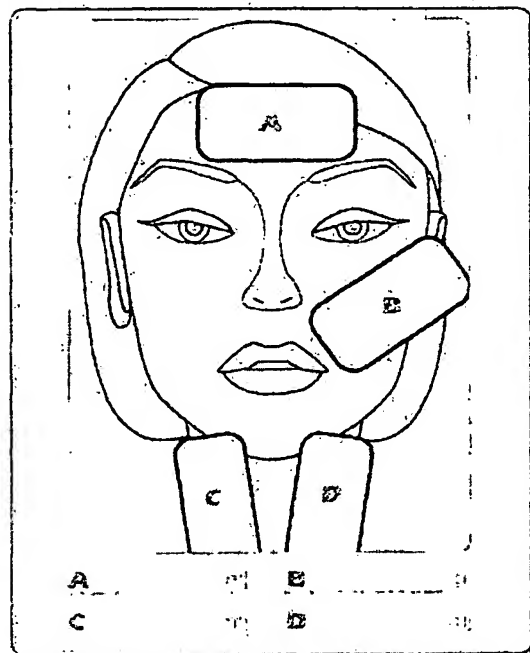

2C

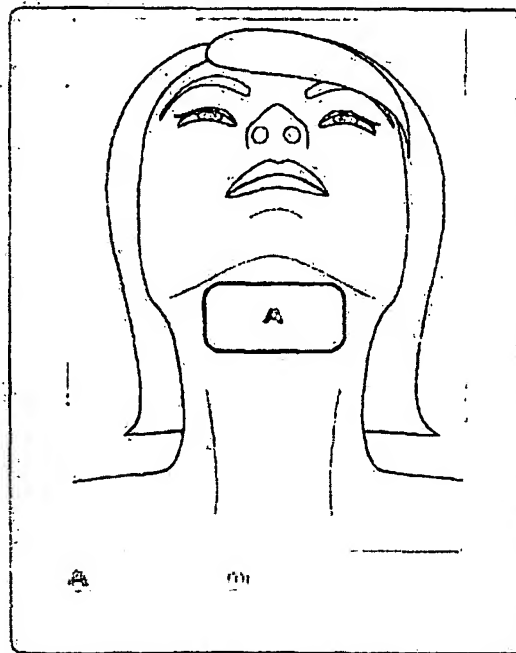

Number of cycles per session: .....

Application time per each electrode: .....



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

---

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA  
ELECTRONIC ENGINEERING SRL**

---

**Patient No. 2**

The picture "2 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...  
during the period between ...10/01/2005..... and .....10/02/2005.....

Electrodes under the following configuration have been used:

- Configuration 1A (electrodes A, B, C, D)

for each electrode it has been set:      density of current ...2.4..... mA / cmq  
   density of power ...100.....mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The pictures "2 After" have been taken at the end of the therapy, and they show considerable improvement with respect to the initial pictures.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine


**Files enclosed:**

- Pictures "2 Before" and "2 After"
- Drawings of the electrodes configurations as used

TEXON AGE


2

Before

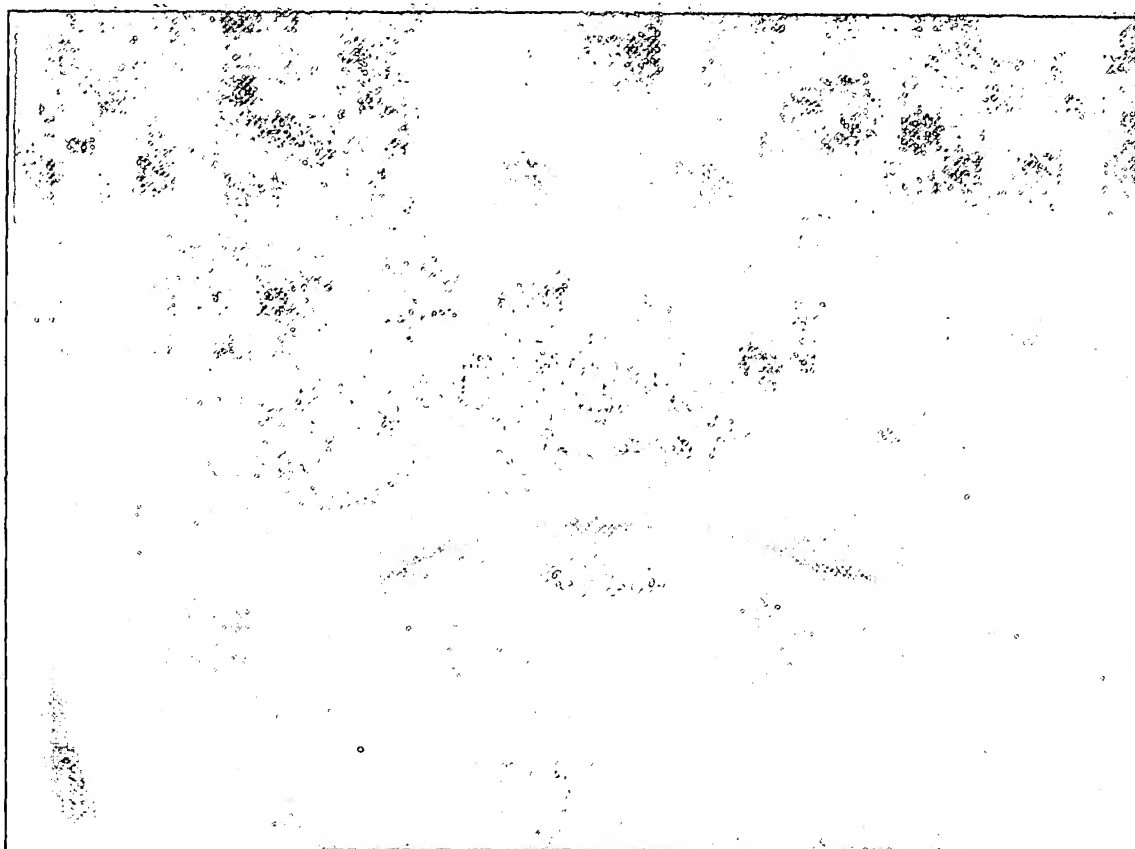


After


# Electrodes setting schemes ELECTRODE nr. 1 (Eye Contour)

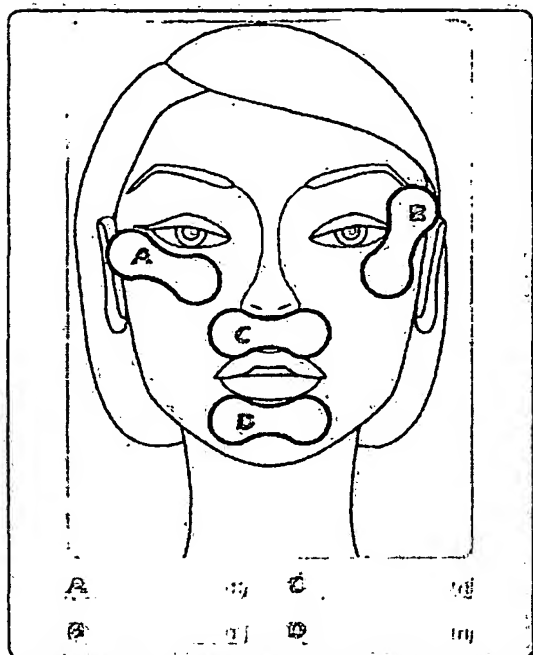

1A

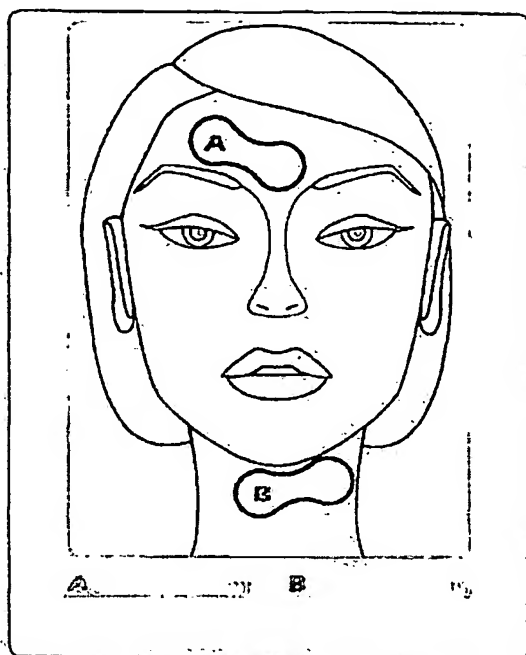

1B

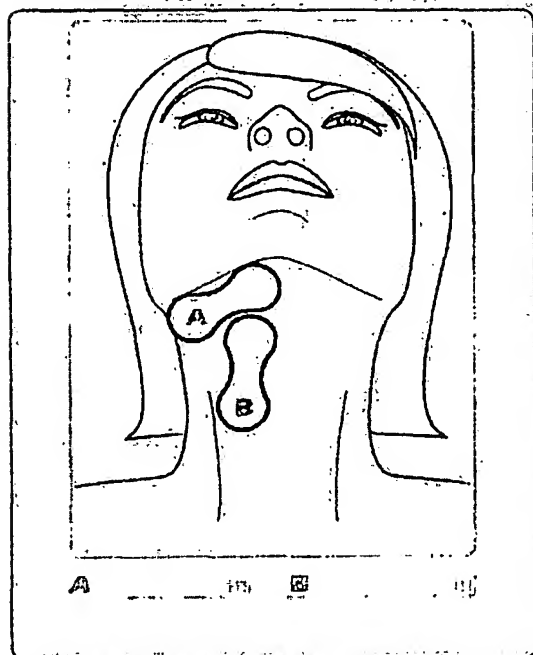

1C

Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

---

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA  
ELECTRONIC ENGINEERING SRL**

---

**Patient No. 3**

The picture "3 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ... 05400106...  
during the period between ... 10/01/2005..... and ..... 10/02/2005.....

Electrodes under the following configuration have been used:

- Configuration 1A (electrodes A, B)

for each electrode it has been set:      density of current ... 2.4 ..... mA / cmq  
   density of power ... 100 ..... mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "3 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

  
.....

**Files enclosed:**

- Pictures "3 Before" and "3 After"
- Drawings of the electrodes configurations as used

TELEA ELECTRONIC ENGINEERING S.r.l.  
Via Leopardi, 23 – 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa – VI-  
Uff. Commerciale: +390444239519 – +390444918179 – Uff. Amministrazione: +390444919322  
FAX: +390444914319 E-mail: [info@vesalius.it](mailto:info@vesalius.it) – <http://www.vesalius.it>  
Codice Fiscale e Partita I.V.A. 01969870243 – Cap. Soc. Euro: 63.191,08 I.V. – Reg. Impr. N. 21604




FEYNMAN 3GE

3


Before

After


Electrodes setting schemes  
**ELECTRODE nr. 1 (Eye Contour)**

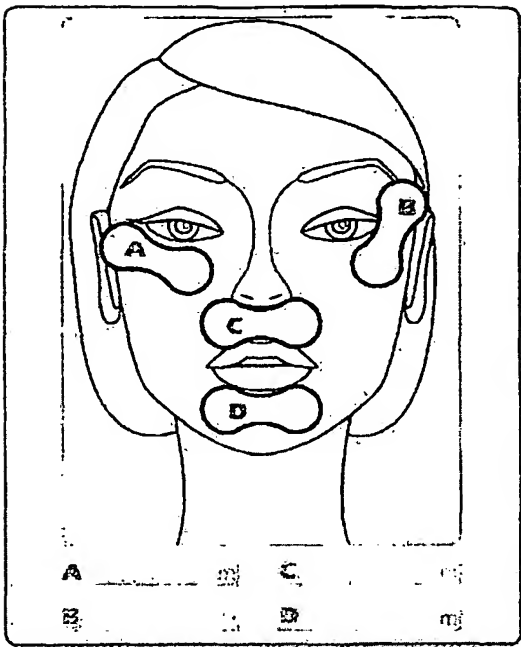

1A

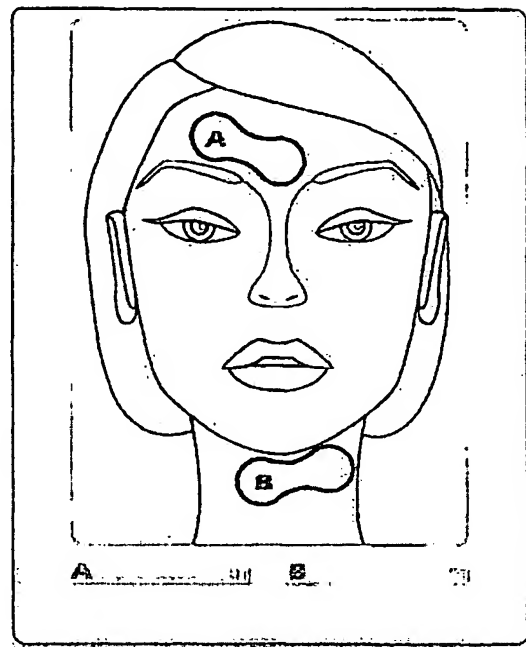

1B

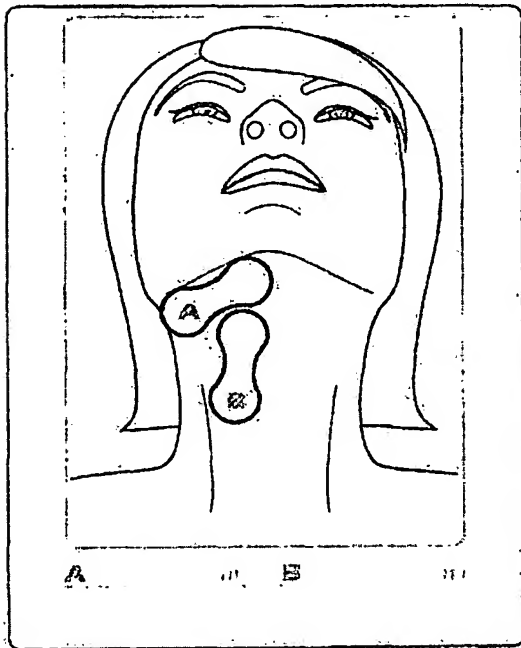

1C

Number of cycles per session: .....

Application time per each electrode: .....



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL**

**Patient No. 4**

The picture "4 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106...  
during the period between ...10/01/2005... and ...10/02/2005...

Electrodes under the following configuration have been used:

**- Configuration 3A (electrode C)**

the electrode has been set: density of current ...2.4... mA / cm<sup>2</sup>  
density of power ...100... mWatt / cm<sup>2</sup>

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "4 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

**Files enclosed:**


- Pictures "4 Before" and "4 After"
- Drawings of the electrodes configuration as used

TEXON 362

4

Before

After

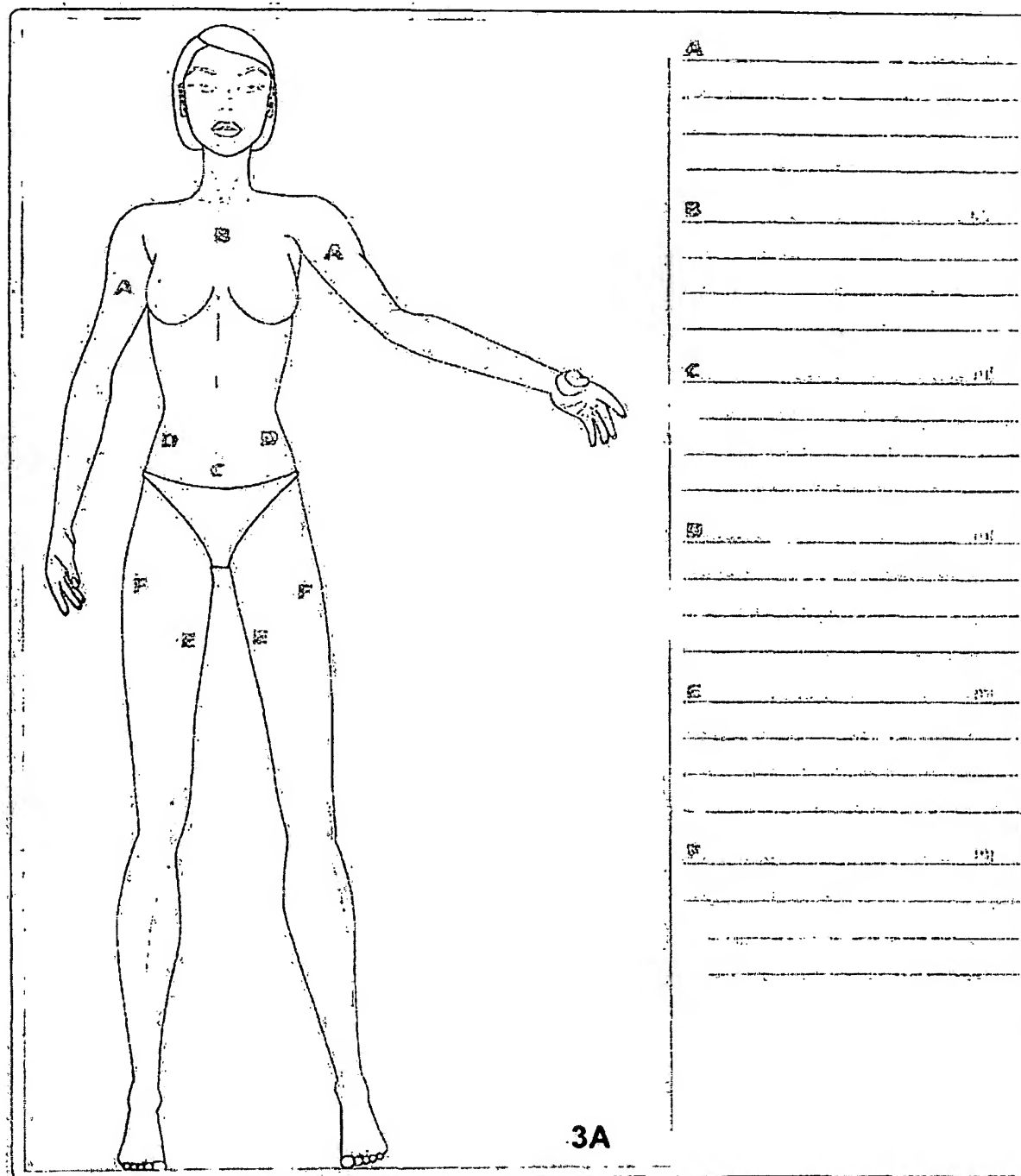

Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

## Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

### Patient No. 5

The picture "5 Before" has been taken before the Rexion-age therapy.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106... during the period between ...15/02/2005..... and .....15/03/2005.....

Electrodes under the following configuration have been used:

- Configuration 3A (electrode C)

the electrode has been set: density of current ...2.4..... mA / cmq  
density of power ...100.....mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "5 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the Rexion-age therapy.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

**Files enclosed:**

- Pictures "5 Before" and "5 After"
- Drawings of the electrodes configuration as used


TELEA ELECTRONIC ENGINEERING S.r.l.  
Via Leopardi, 23 – 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa – VI  
Uff. Commerciale: +390444239519 – +390444918179 – Uff. Amministrazione: +390444919322  
FAX: +390444914319 E-mail: [info@vesalius.it](mailto:info@vesalius.it) – <http://www.vesalius.it>  
Codice Fiscale e Partita I.V.A. 01969870243 – Cap. Soc. Euro 63.191,08 I.V. – Reg. Impr. N. 21604



Before



After

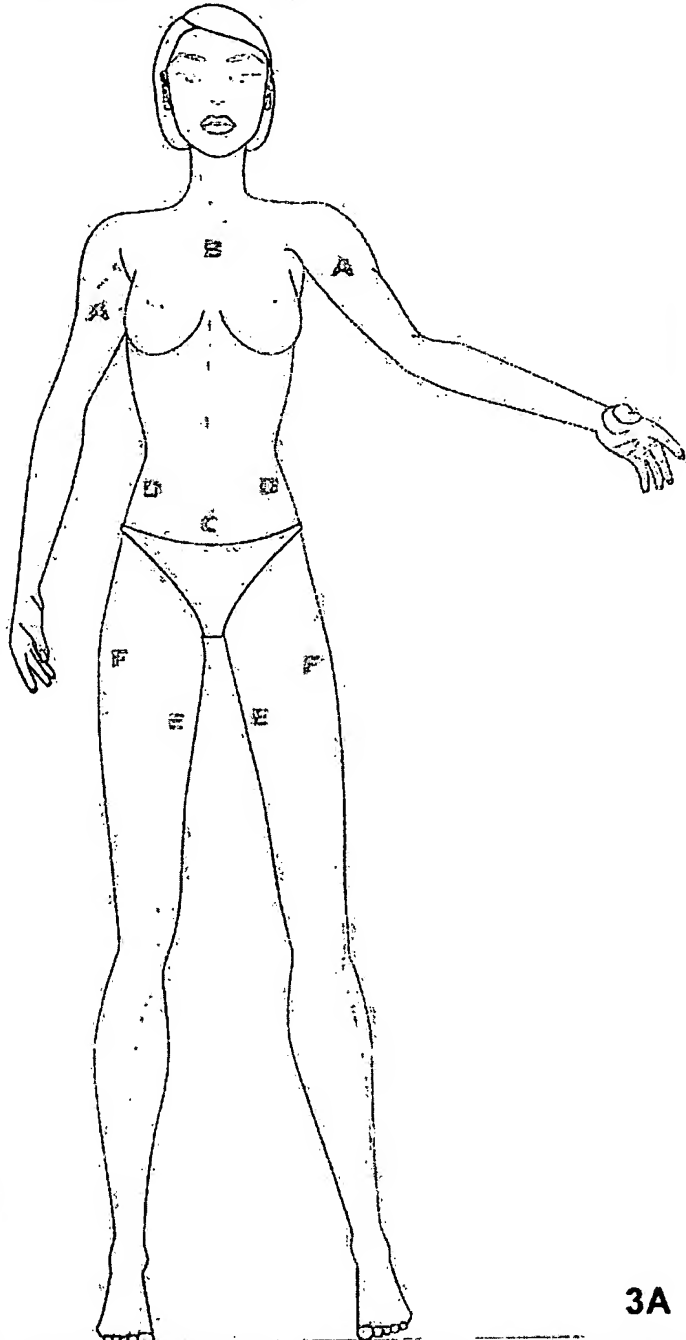

A	_____
B	_____
C	_____
D	_____
E	_____
F	_____
G	_____
H	_____

**3A**

Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_





ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI-EN 46001

## Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL

### Patient No. 6

The picture "6 Before" has been taken before the REXON-AGE therapy.

The therapy has been made with the REXON-AGE unit, whose number plate code is: ...05400106... during the period between ...15/02/2005... and ...15/03/2005...

Electrodes under the following configuration have been used:

- Configuration 3A (electrode C)

the electrode has been set: density of current ...2.4... mA / cmq  
density of power ...100... mWatt / cmq

For the configuration, during each session, every electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.


The picture "6 After" has been taken at the end of the therapy, and it shows considerable improvements with respect to the initial picture.

The improvement obtained has been constant during the 12 months after the end of the REXON-AGE therapy.


Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

### Files enclosed:

- Pictures "6 Before" and "6 After"
- Drawings of the electrodes configuration as used

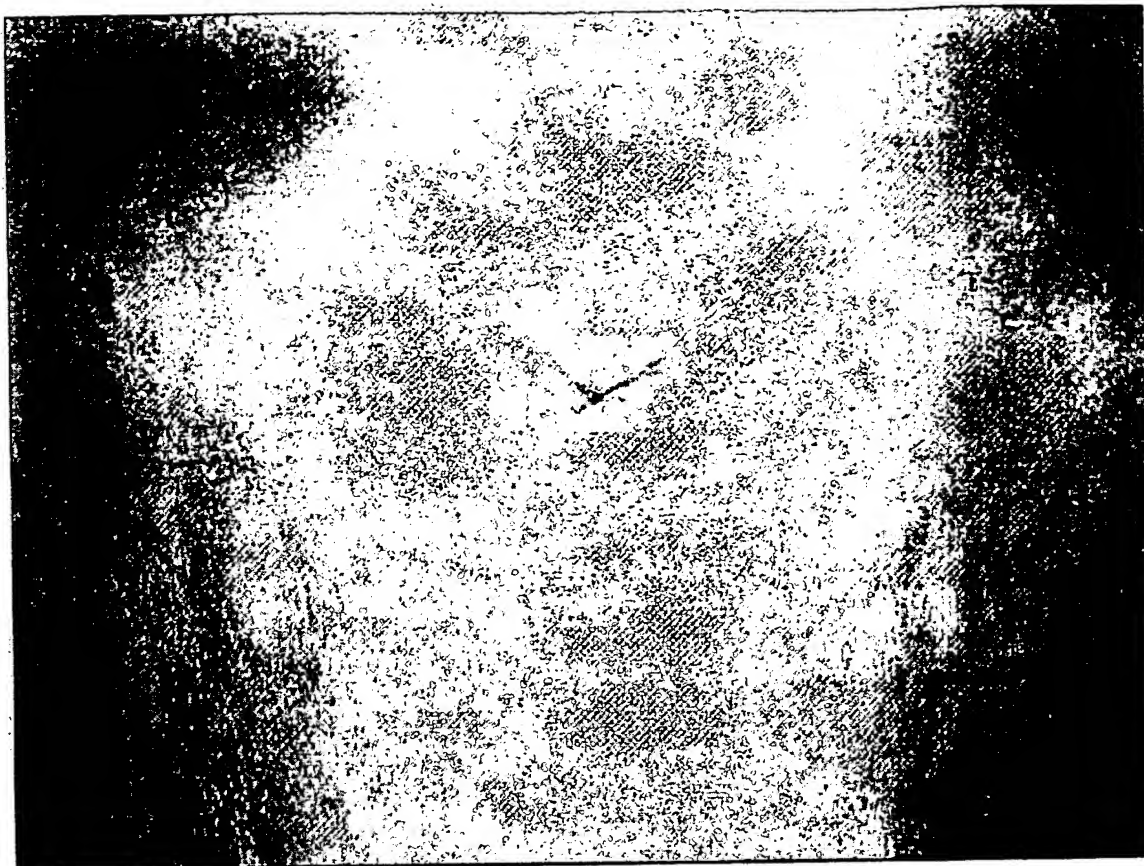


After

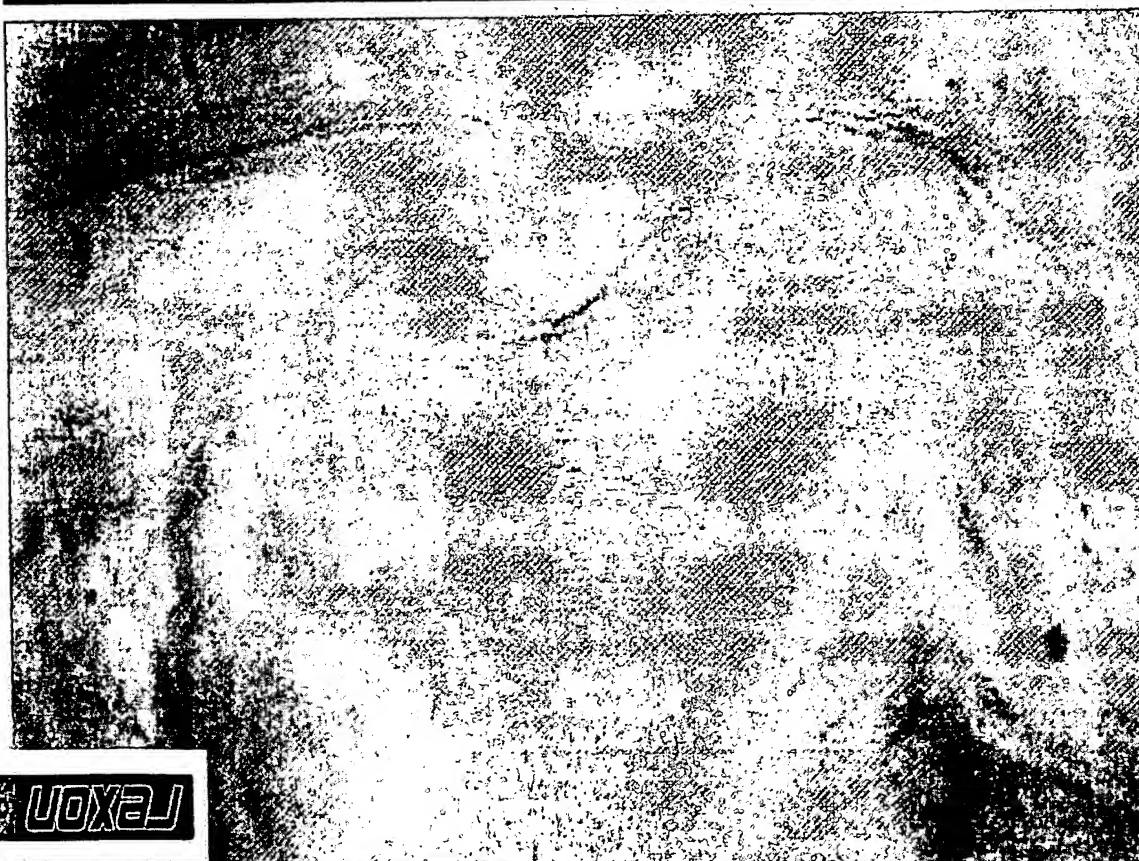

Before

9

Lexon 262

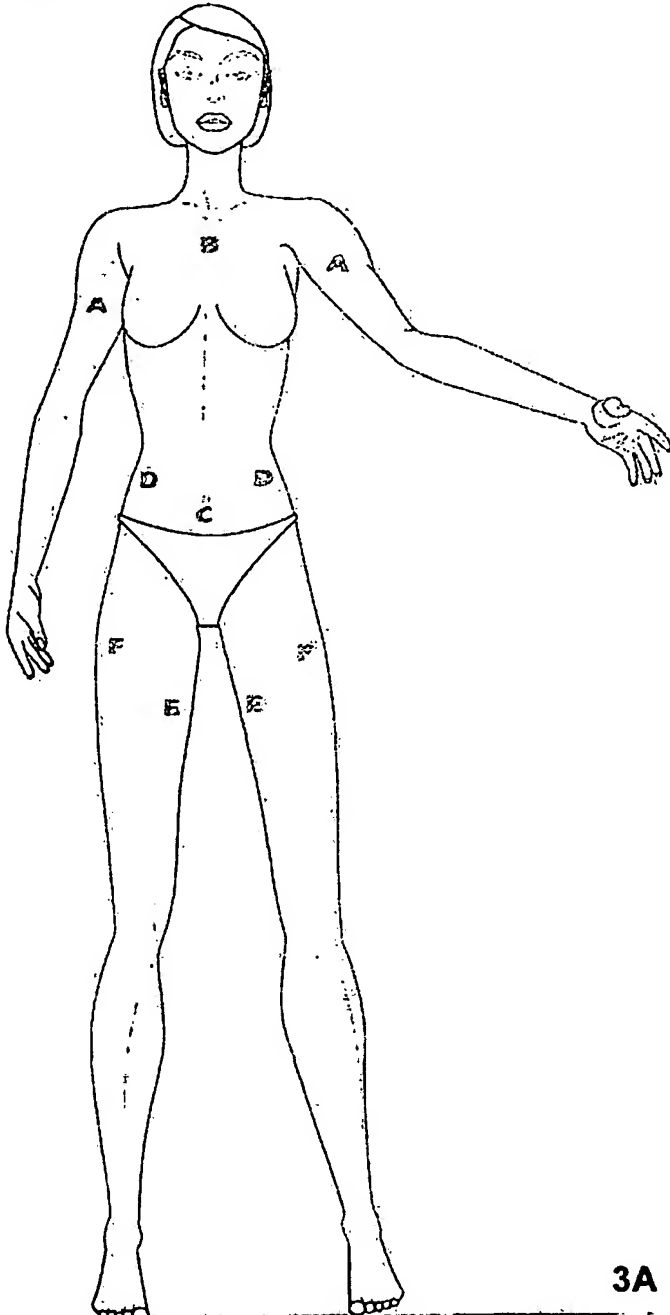

**3A**

<b>A</b>	_____
_____	_____
_____	_____
<b>B</b>	_____
_____	_____
_____	_____
<b>C</b>	_____
_____	_____
_____	_____
<b>D</b>	_____
_____	_____
_____	_____
<b>E</b>	_____
_____	_____
_____	_____
<b>F</b>	_____
_____	_____
_____	_____

Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_



ELETTRONICA  
BIOMEDICA  
E INDUSTRIALE

UNI EN ISO 9001:2000 UNI CEI EN 46001

---

**Clinical tests, concerning the therapy with the Medical Device REXON-AGE produced by TELEA ELECTRONIC ENGINEERING SRL**

---

**Patient No. 7**

The picture "7 Before" shows a section of a portion of skin and lower tissue removed from the inner part of the left arm of the patient No. 7, to whom the electrode A, according to the configuration 3 A, has then been used.

The therapy has been made with the Rexion-age unit, whose number plate code is: ...05400106... during the period between ... 15/02/2005 ..... and ..... 15/03/2005 .....

For the electrode A it has been set: density of current ...2.4..... mA / cmq  
density of power ...100.....mWatt./ cmq

The electrode has been activated for 4 minutes.

In total 6 sessions have been made 3 days apart one another.

The picture "7 After" shows a section of a portion of skin and lower tissue removed from the inner part of the left arm of patient No. 7 at the end of the therapy.


From the picture "7 After" we can see that there has been an important improvement in the quantity of collagen and elastine, due to an increased stimulation on the fibroblasts from the energy transmitted. As a result of the therapy the patient's skin obtained a consistent tightening up and a visible reduction of wrinkles.

Dr. Claudio Trainotti  
Doctor and Lecturer in Aesthetic Medicine

**Files enclosed:**


- Pictures "7 Before" and "7 After"
- Drawings of the electrodes configuration as used

TELEA ELECTRONIC ENGINEERING S.r.l.  
Via Leopardi, 23 - 36050 QUINTO VICENTINO Loc. Villaggio M.te Grappa - VI-  
Uff. Commerciale: +390444239519 - +390444918179 - Uff. Amministrazione: +390444919322  
FAX: +390444914319 E-mail: [info@vesalius.it](mailto:info@vesalius.it) - <http://www.vesalius.it>  
Codice Fiscale e Partita I.V.A. 01969870243 - Cap. Soc. Euro 63,191,08 I.V. - Reg. Impr. N. 21604

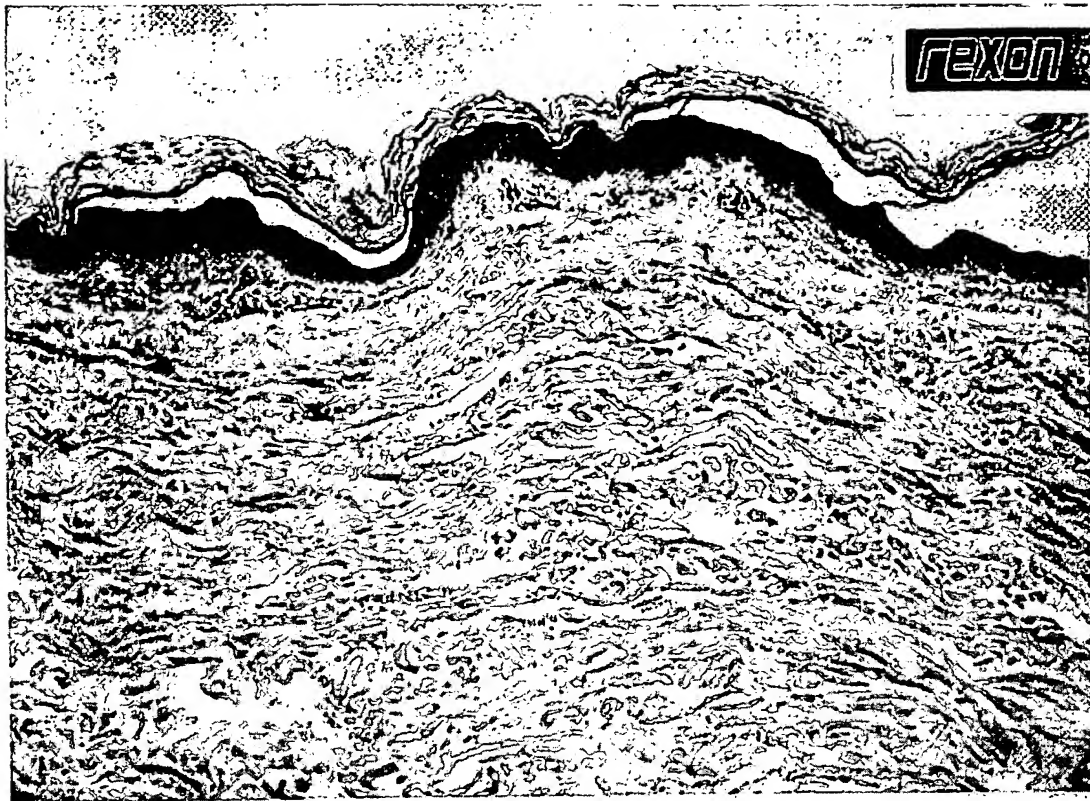


7

Before

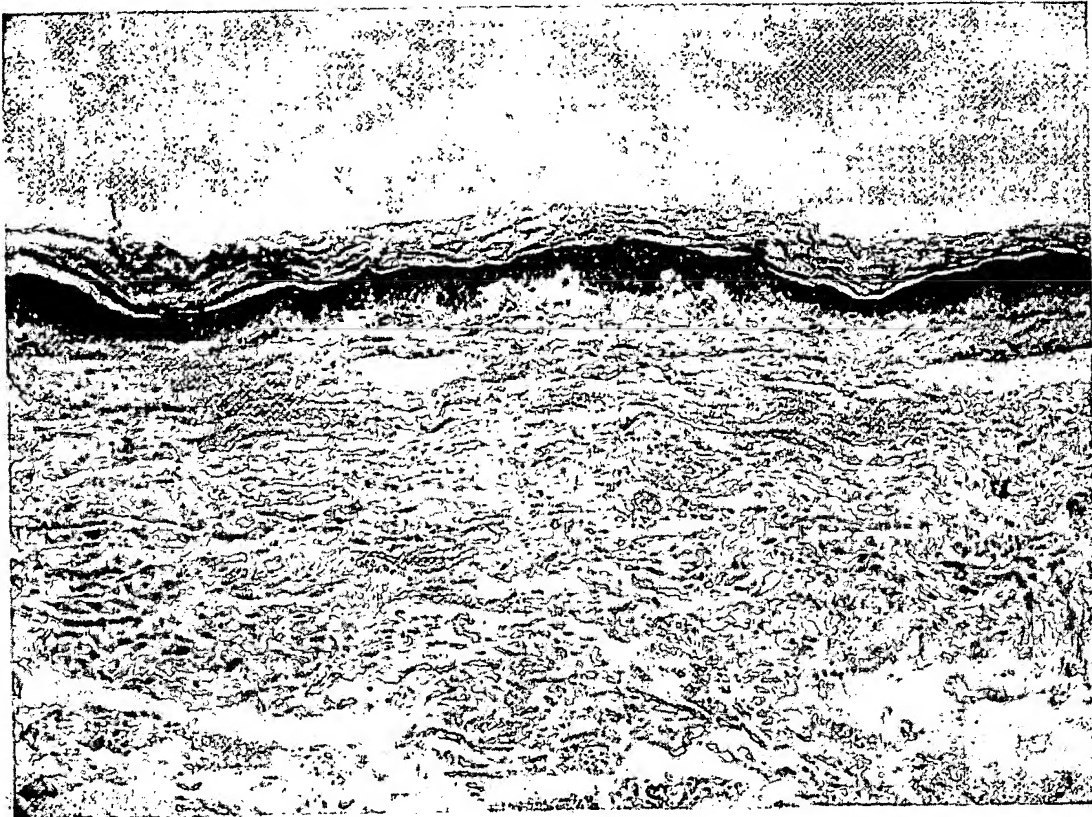

After

Electrodes setting schemes  
**ELECTRODE nr. 3 (Body-front)**

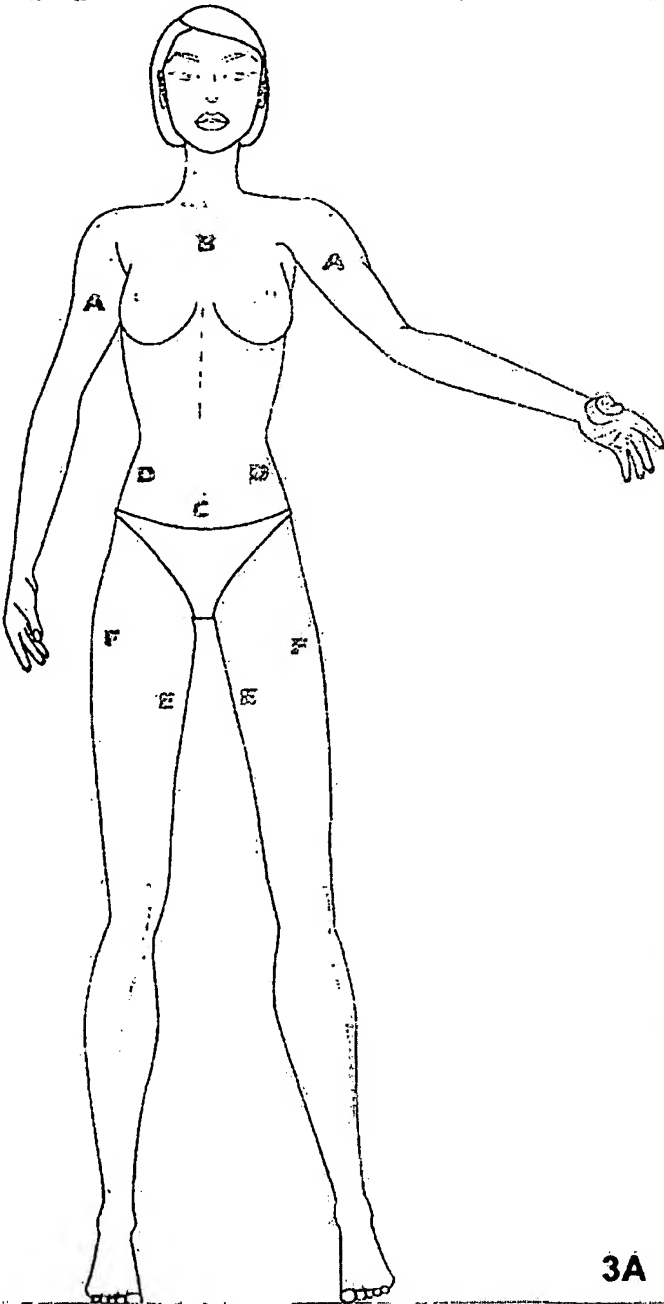	<b>A</b>	_____
	_____	_____
	_____	_____
	<b>B</b>	_____
	_____	_____
	_____	_____
	<b>C</b>	_____
	_____	_____
	_____	_____
	<b>D</b>	_____
	_____	_____
	_____	_____
	<b>E</b>	_____
	_____	_____
	_____	_____
	<b>F</b>	_____
	_____	_____
	_____	_____
	<b>3A</b>	

Number of cycles per session: \_\_\_\_\_

Application time per each electrode: \_\_\_\_\_